


## **Installing Abstract Plus Version 3.1.1.7 Production Release**

To install the Abstract Plus version 3.1.17 production program, complete the following steps:


1. Double-click **setup.exe**

**Result:** The **Abstract Plus Setup Wizard** opens.

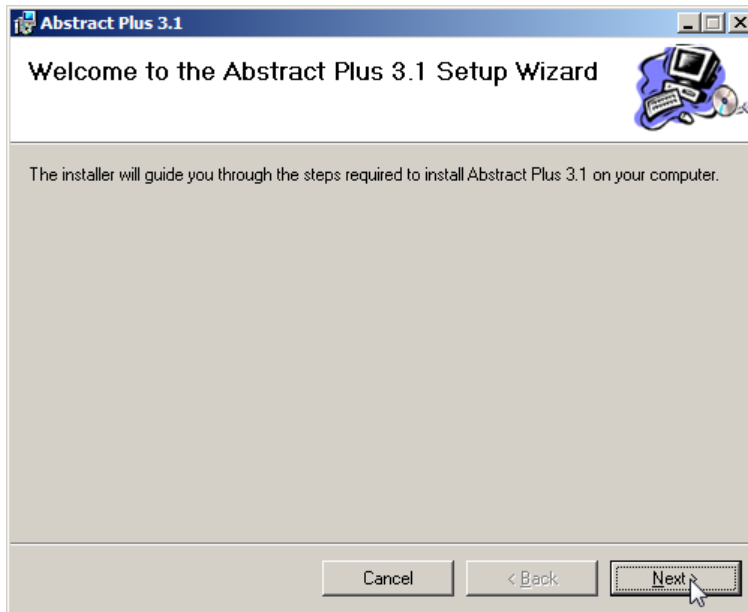


2. Click **Next** to continue installation.

**Result:** The **Select Installation Folder** window opens.

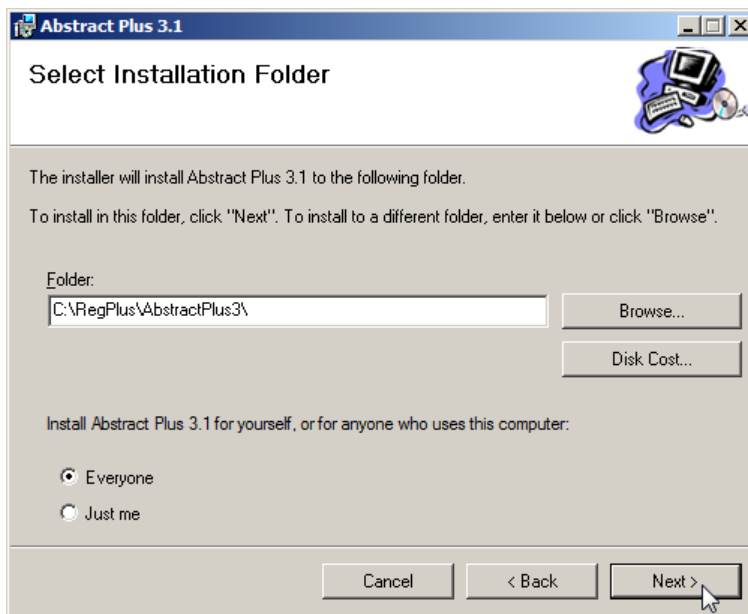


3. Click **Next** to accept the default installation folder.

**Result:** The **Confirm Installation** window opens.

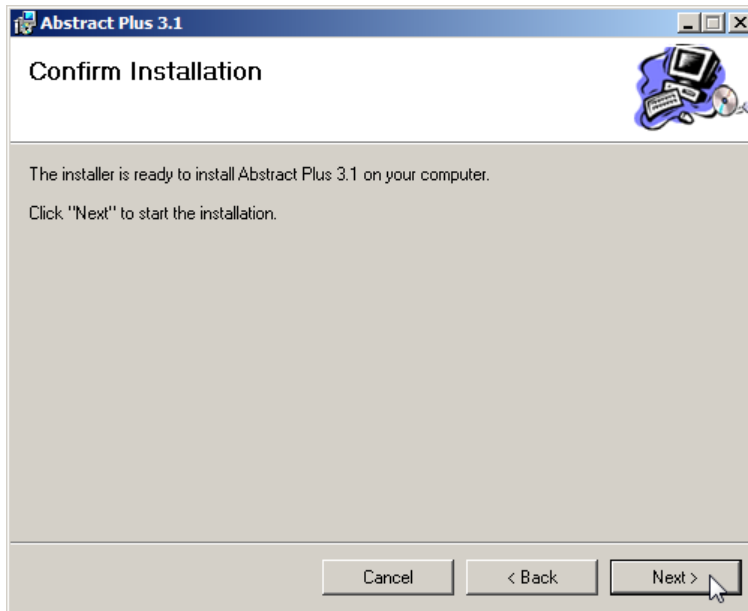


4. Click **Next** to confirm you would like to install Abstract Plus.

**Result:** The installation process begins.

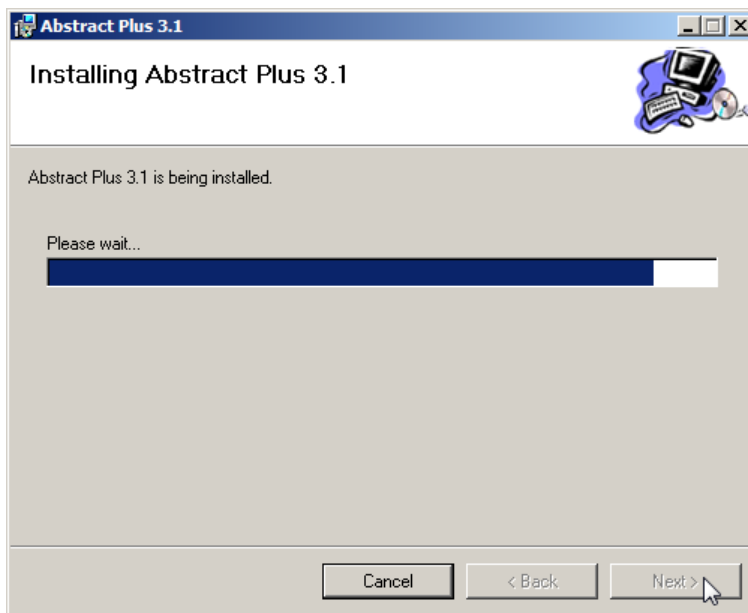

5. When installation has completed, the **Installation Complete** window appears.

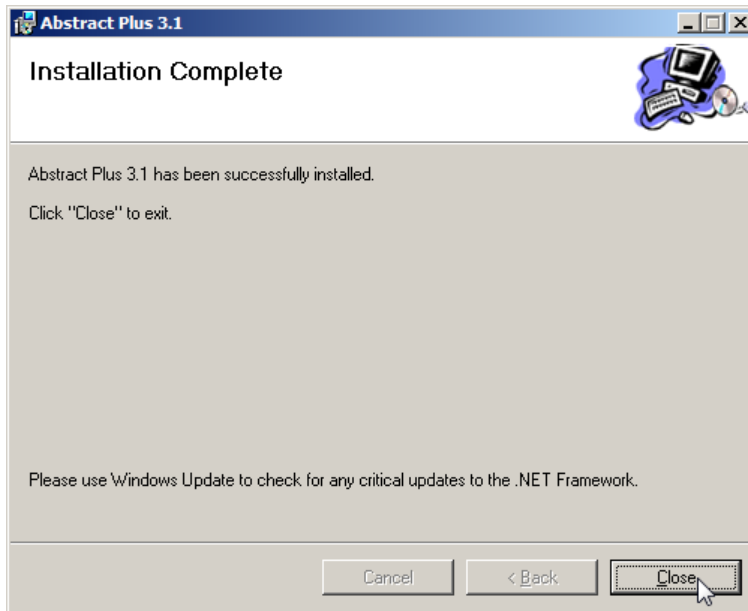


6. Click **Close**.

### **Setting Permissions for User for the Abstract Plus MDBS folder**


**WARNING:** For installations on Vista and Windows 7 operating systems, **Full Control permissions MUST be granted to the user for the Abstract Plus MDBS folder**. If permissions are not set, the database backup and restore features of the program will not work correctly, and may result in data loss in the case of a corrupted database.

To grant Full Control permission to the Abstract Plus MDBS folder please complete the following steps:

1. Launch Windows Explorer (can right-click on the Start menu and select Explore).
2. Locate the Abstract Plus MDBS folder, in **C:\RegPlus\AbstractPlus3\MDBS**
3. **Right-click the MDBS folder**, and select **Properties**.

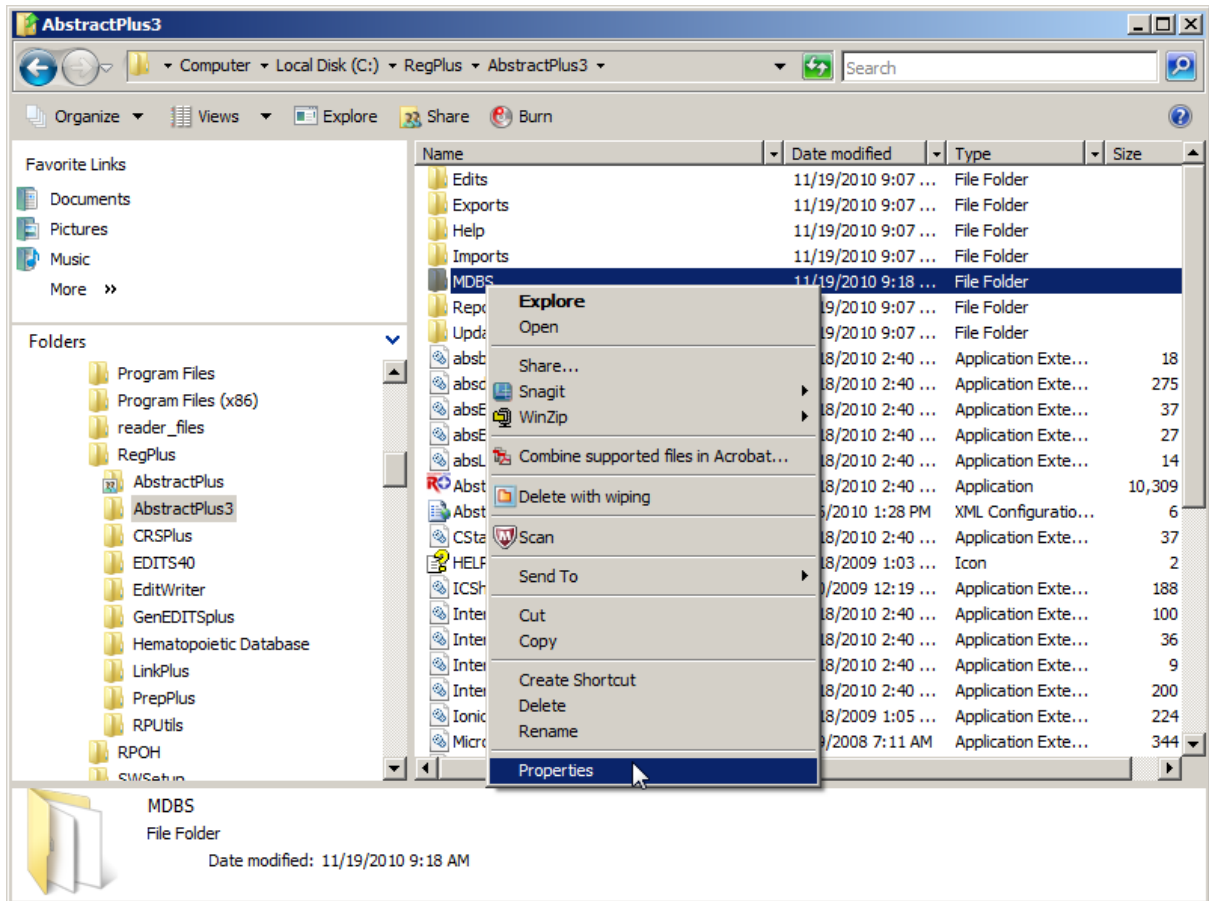


**Result:** The **MDBS Properties** window opens.

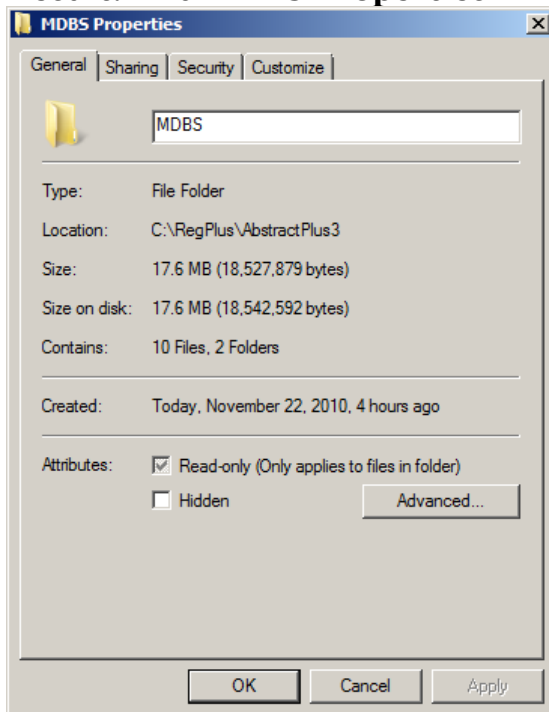


4. Click the **Security** tab.

**Result:** The **Security** tab window opens.

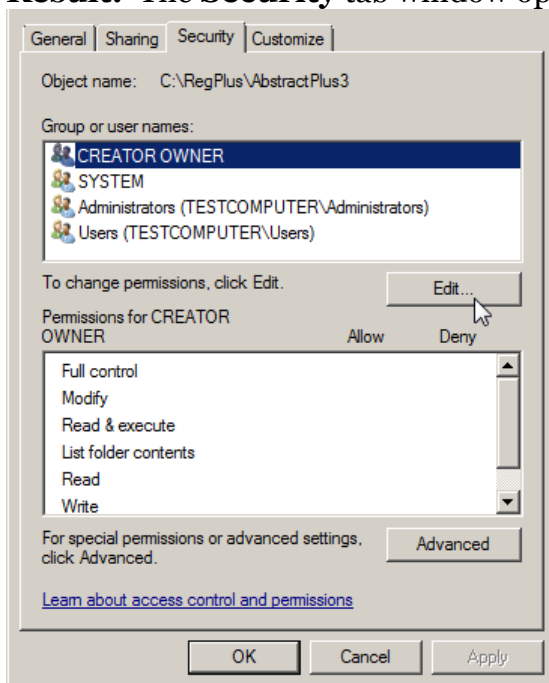

5. Under Group or user names, Click **Edit**.


**Result:** The **Security** tab permissions window opens.

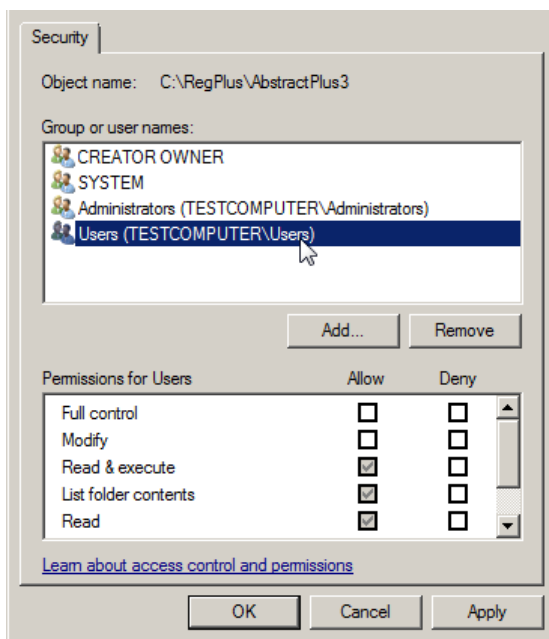


6. Select the Users account to highlight it.

**Result:** The permissions for the Users account is displayed in the lower window


7. In the lower, **Permissions for Users** window, check the Allow check box to the right of the Full control permission (if full control permissions have already been granted, no further action is required).

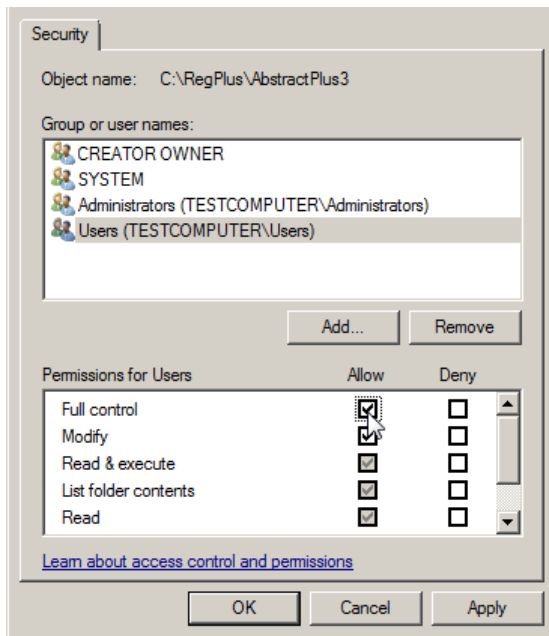

- In the lower, **Permissions for Users** window, check the Allow check box to the right of the Full control permission (if full control permissions have already been granted, no further action is required).

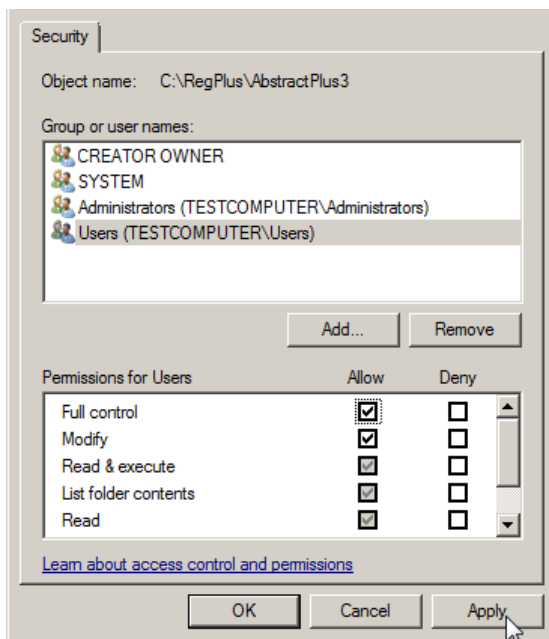


- Click **Apply**.  
**Result:** Full permissions are granted and displayed for the Users account
- Click **OK** to close the **Security** tab permissions window.

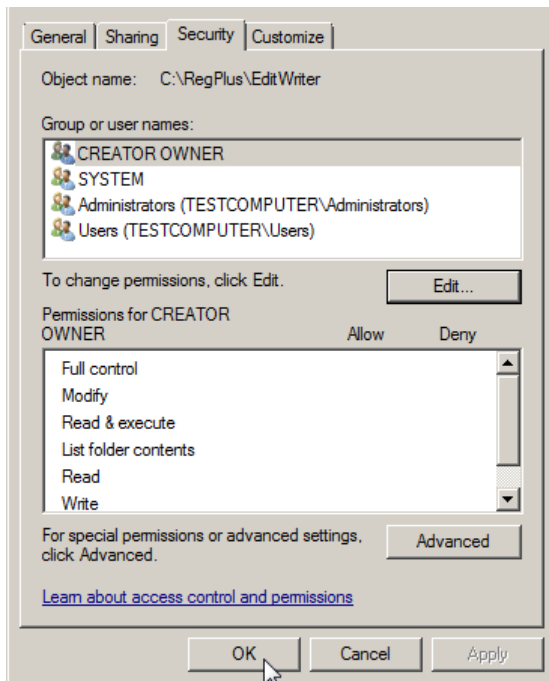

**Result:** You are **returned to the Security tab** window.

11. Click **OK** to close the Security tab window.

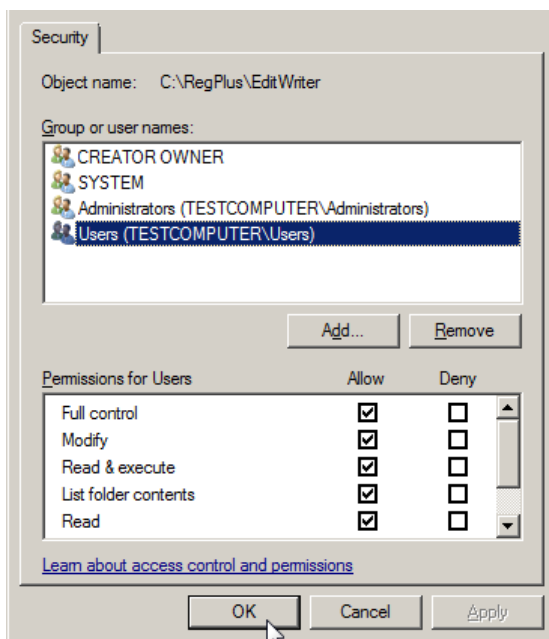

**Result:** You are **returned to Windows Explorer** with the full permissions set.